

CAWO CR Cassettes



user-friendly and durable CR Cassette



Solutions

The CAWO CR Cassette contains the photostimulable phosphor imaging plate.

The CAWO CR Cassette and plate are compatible with existing X-ray tables.

The exposure equipment and routines do not have to be modified when switching from conventional to digital imaging.



User friendly

The CAWO CR Cassette is lightweight and very comfortable to use. It is made of synthetic material and it is warm and pleasant to the touch, minimizing patient discomfort.

Practical maintenance

The life of CAWO CR cassettes is further enhanced by a design which is easily repairable. The closing mechanism and the hinge are easily replaced with the help of a simple screwdriver and are available from CAWO as replacement parts.




Designed to last

The synthetic material of the CAWO CR Cassette has a special ribbed structure that provides maximum rigidity. This makes the CAWO CR Cassette almost indestructible. The corners are rubber protected against accidental dropping.




Antistatic inner lining

The inner lining of the CAWO CR Cassette body is made of felt. This ensures a high degree of protection against electrostatic charging, dust collection and mechanical damage on the imaging plate.

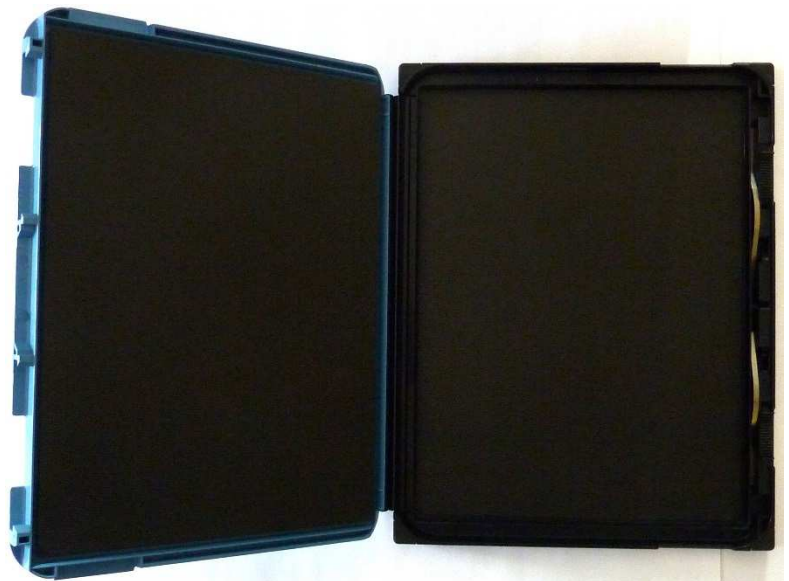

Quality in all situations

Novodur® is a particularly temperature-resistant material, bearing up well under extremes in climatic conditions. Another plus for the reliability and the life expectancy of the CAWO CR cassette. This applies in whatever part of the world you want to use them and in whatever situations.

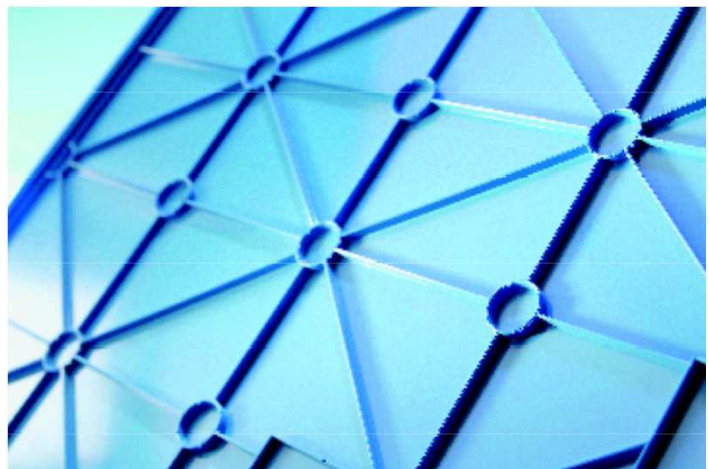

Metric

13 x 18 cm	24 x 30 cm
15 x 30 cm	18 x 24 cm
30 x 40 cm	18 x 43 cm
35 x 35 cm	20 x 40 cm
35 x 43 cm	

Imperial

8 x 10 in	14 x 14 in
10 x 12 in	14 x 17 in

Weight (35 x 43 cm):

1.6 kg

CAWO CR cassettes comply with the requirements of the Council Directive 93/42/EEC (Medical Devices) and carry the CE mark. They are in conformity to the international standard DIN ISO 4090: 2004-08.



CAWO Version 2019-02

A wide range of secondary radiation grids, aluminium **CAWO MMS-FA Cassettes** for long formats and grid applications, blue and green-emitting intensifying screens, accessories, equipment ensure tailor made solutions for many radiologic applications.

Business Contact:

CAWO Solutions

P.O. Box 11 29

D - 86521 Schrobenhausen

Phone: (++49) 8252 9109-21

Internet: www.cawo.com

E-mail: cawo@agfa.com

These products are designed for medical applications only. They shall not be used for Non-Destructive-Testing (NDT) applications.

CAWO Solutions **CE**

Agfa NV - Septestraat 27 - 2640 Mortsel - Belgium