

CAWO Screen Cleaner



Solutions

CAWO Screen Cleaner

An easy-to-use, isopropyl-based cleaner that ensures longer life and higher performance for intensifying screens.

TECHNICAL SPECIFICATIONS

Solvent

- Isopropyl, contains an Anti static agent

Contains

- 4 x 500 ml

Safety

- Highly flammable
- Keep container closed
- Keep away from sources of ignition


Safety datasheet available on request

Protect your investment

Dust and dirt build up on intensifying screens, leading to image artefacts. To avoid this, screens need cleaning at regular intervals. Using this cleaner, specially designed for intensifying screens, will help ensure them the longest service life.

Simple to use

CAWO Screen Cleaner is simple to use: moisten a piece of gauze with the cleaner, wipe gently across the screens surface and wait for any remaining solvent to evaporate.

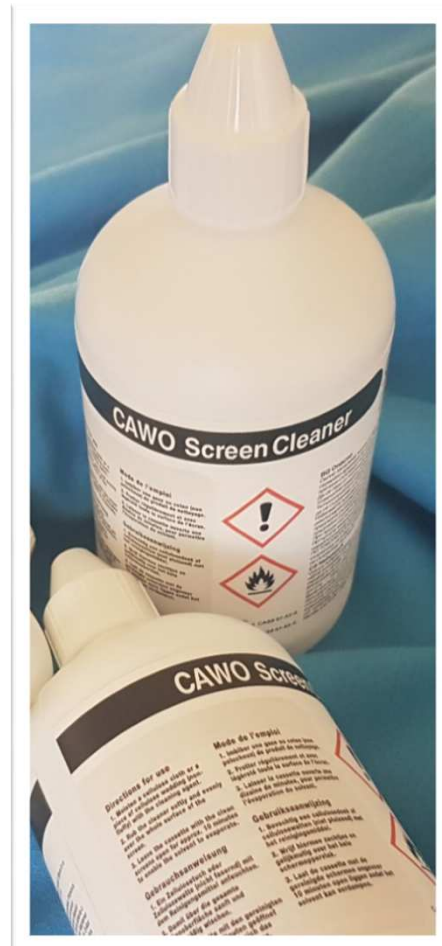

Business Contact:
CAWO Solutions
P.O. Box 11 29
D - 86521 Schrobenhausen
Phone: (+49) 8252 9109-21
Internet: www.cawo.com
E-mail: cawo@agfa.com



Solutions

Agfa NV - Septestraat 27 - 2640 Mortsel - Belgium